



# Innovative. Smart. Trusted.

We're your partner for sterilization products.



Monitoring Products for Vaporized Hydrogen Peroxide

VH202

Select a system that works in all low temperature sterilization systems. Cintel Medical products for monitoring Vaporized Hydrogen Peroxide (VH<sub>2</sub>O<sub>2</sub>) sterilization processes provide fast, visible, consistent and reliable results offering peace of mind for assurance of effective sterilization.

## SporView® Self-Contained Biological Indicators

- Compatible with all low temperature sterilization systems, including Sterrad®
- *Geobacillus stearothermophilus* product configurations available
- Exclusively formulated Trypticase® Soy Broth (TSB) modified with pH indicator
- Visually distinct positive results from purple to yellow
- Chemical Indicator on label provides immediate indication of exposure
- Smaller size, 8.4 mm x 45.5 mm, designed to fit in challenging sterilization areas
- Reduces the possibility of contamination and minimizes the potential for false positives
- Incubation time of 24 hours
- Compliant with ISO 11138-1





## Dry Block Incubators

- Designed for convenient and reliable incubation of SporView® SCBIs
- Electronically-heated aluminum block
- Holds up to 10 or 15 units
- Adjustable temperature (5°C to 70°C) and incubation time (0 to 99h 59m)
- Reaches set point temperature in 15 minutes or less
- Transparent lid included to ensure consistent temperature and to protect the units

## Chemical Indicator Labels

- Provide a fast and reliable visual check on items that were exposed to VH<sub>2</sub>O<sub>2</sub>
- Consistent color transition from red to blue
- Permanent pressure-sensitive adhesive
- Room temperature storage (15°C to 30°C)
- Meet performance specifications as outlined in ISO 11140-1



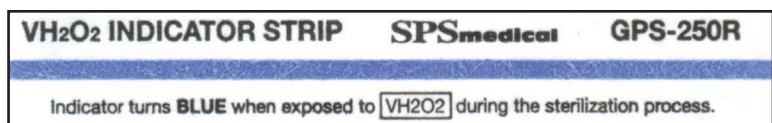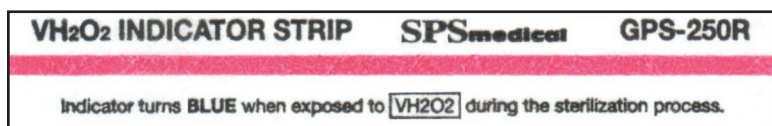

## Chemical Indicator Strips

- Certified to meet chemical indicator specifications for pack monitoring when used in VH<sub>2</sub>O<sub>2</sub> sterilization processes
- Provides a distinct color change from red to blue upon exposure to VH<sub>2</sub>O<sub>2</sub> processes
- Demonstrates exceptional performance in the most widely used VH<sub>2</sub>O<sub>2</sub> sterilizers
- Manufactured with synthetic material
- ISO 11140 - 1 Type 4 Compliant


## Indicator Tape

- Designed to aid in distinguishing processed from unprocessed loads
- Consists of a strong custom-made white polyester substrate with a pressure sensitive adhesive which adheres to Tyvek®, plastic, woven, treated woven and non-woven wraps
- Easy to remove adhesive allows for secure seals
- Manufactured free of latex & lead
- Available color options from red to yellow or red to blue when exposed to 2.3 mg/L ± 0.4 mg/L of VH<sub>2</sub>O<sub>2</sub> at 50°C for 6 minutes exposure time or longer
- Each shipment includes a Certificate of Analysis for assurance of product quality and performance
- ISO 11140-1 Type 1 Compliant



## Tyvek® Sterilization Pouches

- Single-use pouches come in self-seal, heat-seal and tubing for custom length
- Manufactured with Tyvek® material for superior tear strength and puncture resistance
- Includes a chemical indicator for verification of exposure
- Maintains instrument sterilization for up to 24 months
- Pouches and tubings feature a convenient built-in green indicator that will turn yellow when exposed
- Meets ISO 11607 standards for pouch performance

PRODUCT IMAGE	DESCRIPTION	CATALOG NO.*	UOM	SOLD PER
	<b>SporView® VH<sub>2</sub>O<sub>2</sub> Self-Contained Biological Indicator</b>	VH50	50 Vials/Box	Box
	<b>Dry Block Incubator</b>	NDB-060	1 Ea.	Ea.
	<b>Chemical Indicator Label</b>	GPL-2000R	2000/Roll	Roll
	<b>Chemical Indicator Strip</b>	GPS-250R	250/Box	Box
	<b>Universal Lot Control Label</b>	URL-Series	1125/Roll, 12 Rolls/Box	Box
	<b>Chemical Indicator Tape</b>	NMTAP-H	48 Rolls/Case	Roll
	<b>Sterilization Pouch</b>	TSP-Series THP-Series	200/Box 4-10 Boxes/Case	Box
	<b>Sterilization Tubing</b>	TST-Series	100ft/Roll 10 Rolls/Case	Roll

## Protect with the pros you can trust: Cantel

A leading global company dedicated to delivering innovative infection prevention and control and sterility assurance solutions, Cantel offers a broad line of high-quality Vaporized Hydrogen Peroxide (VH<sub>2</sub>O<sub>2</sub>) products. Each is backed by a team of dedicated experts you can count on to be there to help you with compliance – how, when and where you need them – as you work to improve outcomes, enhance safety – and save lives.

\* Additional products available



Unless otherwise specified, all trademarks are property of Crosstex International, Inc. or SPSmedical Supply Corp. Trypticase® is a registered trademark of Becton Dickinson and Company. Tyvek® is a registered trademark of E.I. DuPont de Nemours and Company.

Cantel (Canada)  
90 Gough Road, Unit 1  
Markham, ON L3R 5V5

Toll Free Phone: +1.844.348.5637  
Toll Free Tax: +1.844.348.5637  
Email: info@cantelmedical.ca

[www.cantelcanada.com](http://www.cantelcanada.com)