

ADVANTAGE PLUS™

Endoscope Reprocessing System | **Reprocess**





THE COMPLETE CIRCLE OF PROTECTION

As the global vanguard in infection prevention, only Cantel Medical delivers the Complete Circle of Protection, a full-value, proactive partnership dedicated to helping you remove risk, streamline operational efficiencies and optimize your success.

The Most Advanced Automated Endoscope Reprocessor From the Market Leader

Reprocess a Variety of Endoscopes

- Compatible with all flexible endoscopes
- Validated for reprocessing TEE probes, dilators, cystoscopes and semi-critical rigid scopes

Comprehensive Endoscope Channel Monitoring

- Guarantees proper hookup connection
- Monitors endoscope channels for blockage
- Ensures all channels are fully reprocessed

High Performance Cleaning Cycle

- Automated cycle eliminates variability of manual cleaning
- Optimizes endoscope HLD process
- Documented wash cycle

Automated Leak Testing and Alcohol Flush

- Enables disinfection of an endoscope prior to repair
- Minimizes leak-related repair costs
- Facilitates drying of endoscope channels

Networking and Remote Diagnostics

- Network multiple systems together for easy administration and reporting
- Connect to local or off-site server for centralized data collection and convenient back-ups
- Fast technical support via on-line, remote diagnostics
- Integrates with ENDORA™ Endoscope Track & Trace System

RAPICIDE™ PA High-level Disinfectant

- Proven endoscope compatibility
- Faster endoscope availability
- Reduced worker exposure; environmentally friendly
- Low-temp disinfection and Aldehyde-free



REPROCESS

High-level disinfection is the cornerstone of infection prevention. Reprocessing technologies from Cantel are designed to ensure patient safety by adhering to the strictest of standards for high level disinfection, optimizing workflow efficiency, and delivering versatility by supporting a wide range of scopes.

ADVANTAGE PLUS™

Endoscope Reprocessing System | **Reprocess**

Visual display of real-time reprocessing information



Extensive cycle data management integrates with ENDORA™ Endoscope Track & Trace System

Large, easy loading, independent basins have asynchronous operation and accommodate all flexible endoscopes



High performing cleaning cycle

Remote Diagnostics for efficient servicing



Hands-free lid operation provides ease of use and decreases endoscope cross-contamination risk

Barcode reader for fast and accurate data entry. PC based for scope specific cycle programming



Individual channel monitoring for endoscope blockage and hook-up connection



Convenient transport cassette for endoscope protection and recontamination prevention



Dimensions

59 H x 45 W x 36 D inches (150 H x 114 W x 92 D cm)
 Lid Height: Open - 54" (137 cm), Closed - 38" (97 cm)
 Depth with air compressor: 37.5" (95 cm)

Weight

400 lbs (182 kg)

Electrical Requirements

100-240 VAC ± 10%: single phase, 50/60HZ; 1000 Watts

Electrical Safety Certifications

ETL Intertek Listed #75811
 UL Standard 61010-1
 EN Standard 61010-1
 CSA/CAN Standard C22.2, No. 61010-1

Water Supply

Minimum flow rate: 1.32 GPM (5 LPM) at 35 to 87 PSI (2.4 to 6 bar)
 Water Temperature:
 (35°C ± 2°C) - North America; (25°C ± 2°C) - International
 Water Hardness: less than 12 gpg (200 ppm)
 Pre-filter wall space: 26 H x 38 W x 6 D inch (66 H x 96 W x 15 D cm)

Basin Size


19 W x 17 L x 7 D inch (48 H x 43 L x 18 D cm)

ADVANTAGE PLUS™ Endoscope Reprocessing System Ordering Information - North America

Model Number	Description
ADV-1007	ADVANTAGE PLUS Endoscope Reprocessing System for RAPICIDE™ PA Disinfectant with air compressor
ADV-1008	ADVANTAGE PLUS Endoscope Reprocessing System for RAPICIDE PA Disinfectant without air compressor
CAS-1000	Endoscope Transport Cassette

ADVANTAGE PLUS Endoscope Reprocessing System Ordering Information - International

Model Number	Description
ADV-2001	ADVANTAGE PLUS Endoscope Reprocessing System for RAPICIDE PA Disinfectant with air compressor
ADV-2002	ADVANTAGE PLUS Endoscope Reprocessing System for RAPICIDE PA Disinfectant without air compressor

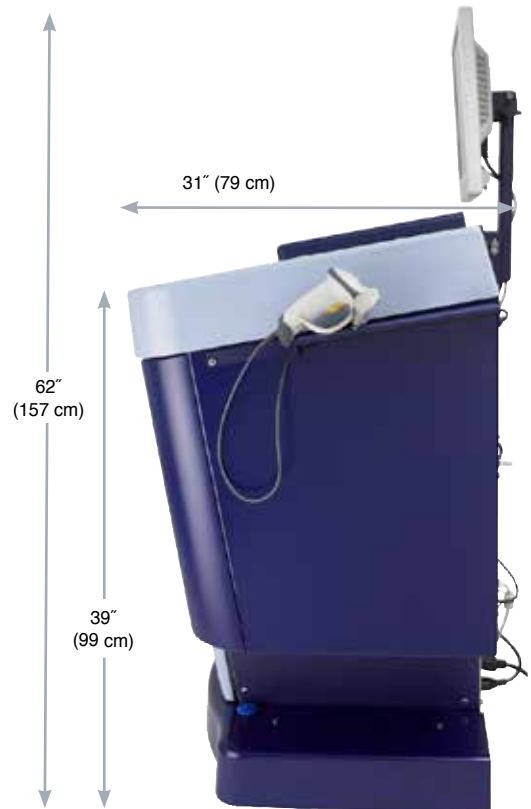

ADVANTAGE PLUS™, ENDORA™ and RAPICIDE™ are trademarks of Medivators Inc.



Cantel Medical (Canada)
 90 Gough Road, Unit 1
 Markham, ON L3R 5V5

Phone: +1.844.348.5636
 Fax: +1.844.348.5637
 Email: info@cantelmedical.ca

www.cantelcanada.com

© 2017 Medivators Inc. P/N 50098-077C-ENCAN