

# VERISCAN™ LT

Automated Endoscope Leak Detection System | **Manual Clean**






## THE COMPLETE CIRCLE OF PROTECTION

As the global vanguard in infection prevention, only Cantel Medical delivers the Complete Circle of Protection, a full-value, proactive partnership dedicated to helping you remove risk, streamline operational efficiencies and optimize your success.

## A Revolution in Endoscope Leak Testing

### Lowers Cross-contamination Risk Due to Undetected Leaks

- Detecting 99% of endoscope leaks<sup>1</sup>
- Utilizes highly sensitive, innovative leak sensing technology
- Unsurpassed automated leak detection, compared to approximately 65% by visual observation methods<sup>2</sup>
- Dry endoscope leak test – no fluid submersion required
- Compatible with all flexible endoscopes
- Detects leaks as small as 0.0008" (~1/3rd the diameter of a human hair)

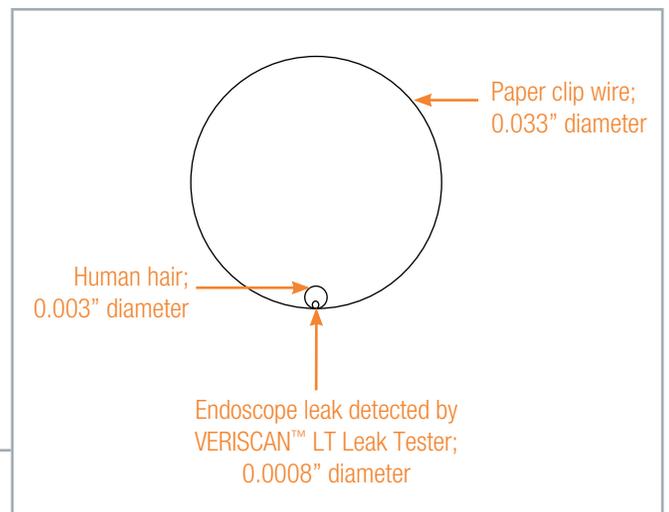

### Ensures Compliance with Endoscope Leak Test Protocols

- Consistent, repeatable endoscope leak testing
- Printable confirmation of leak test results
- Documented endoscope, operator and patient ID information

### Simple to Operate

- Touch Pad control
- Visual display of cycle screen steps
- Quick, simple endoscope connection

*“Detecting a leak early can result in inexpensive repair. It also prevents fluid invasion, which represents more than 50 percent of total repair costs (that can range between \$1,000 - \$8,000), depending on the damage. Identifying and treating scopes is vital to keeping repair costs down, as well as protecting patient health.”*

**Chartier, Keith**  
(2009, April/May). **Keep Your Patients Safe by Checking for Scope Leaks.** *EndoNurse*, V9, 20.

Source:

1. MEDIVATORS VERISCAN LT Endoscope Leak Detection System; Leak Test Accuracy Summary; June, 2011.
2. Verimetrix Company Veriscan Leak Testing Technology Presentation, January 22, 2009.

## MANUAL CLEAN



Manual cleaning is essential in removing bioburden and preparing the endoscope for high-level disinfection. For optimal results, Cantel's manual cleaning products and detergents are designed to be a critical first step to efficiently and effectively eliminating bacteria on endoscopes prior to reprocessing.

# VERISCAN™ LT

Automated Endoscope Leak Detection System | **Manual Clean**

Proven to Drive Down  
Endoscope Repair Costs;  
Detects Endoscope Leaks  
Early When Repairs are  
Minor and Less Costly



Visual display  
steps operator  
through the  
process

Large keypad for  
ease of operation

Simple and quick  
endoscope connection



Documented  
leak test results

VERISCAN LT  
SM 201120001010  
Version: 2.22  
03/30/09 09:36  
Scope S/N: \*\*\*\*\*  
Run ID: 9  
Test Duration: 60 sec  
Test Result: PASS  
Test Setting: STANDARD I  
Knobs Turned  
\*\*\*\*\*



## Available Accessories



Crystal Desiccant



Coalescent Filter



Printer Paper

### VERISCAN™ LT Leak Tester

#### Dimensions

13.125 H x 8.25 W x 19.25 D (inches)  
33.3 H x 21 W x 48.9 D (cm)

#### Unit Weight

15.8 lbs (7.18 kg)

#### Electrical Safety Certifications

UL Standard 61010-1  
CAN/CSA Standard 61010-1

#### Electrical Requirements

100-240VAC, 50/60Hz, .25 amp

#### PSI Range

3.5 – 3.9 PSI (.23 to .26 bar)

#### Warranty

One year

### VERISCAN LT Leak Tester Ordering Information

Part Number	Description
VER-LT	VERISCAN LT Endoscope Leak Tester

### Accessory Ordering Information

Part Number	Description
40100-089	Printer Paper (package of five rolls)
67720-596	Crystal Desiccant (package of one)
MF01-0056	Coalescent Filter (package of one)

### VERISCAN LT Leak Tester Hook-up Ordering Information

Part Number	Description
VER-2501	Install Kit for PENTAX® Endoscopes
VER-2511	Install Kit for KARL STORZ® Endoscopes
VER-2521	Install Kit for FUJIFILM™ Endoscopes
VER-2531	Install Kit for OLYMPUS® Endoscopes



VERISCAN™ is a trademark of MEDIVATORS Inc.  
PENTAX® is a registered trademark of Hoya Corporation  
KARL STORZ® is a registered trademark of Karl Storz & Co.

FUJIFILM™ is a trademark of FUJIFILM Corporation  
OLYMPUS® is a registered trademark of Olympus Corporation

[www.cantelcanada.com](http://www.cantelcanada.com)

#### TO PLACE AN ORDER

p: 1.844.348.5636 f: 1.844.348.5637 e: [orders@cantelmedical.ca](mailto:orders@cantelmedical.ca)

© 2017 Medivators Inc. P/N 50098-032B-ENCAN

